


Chair mounted Holder type



Chair mounted Continental type





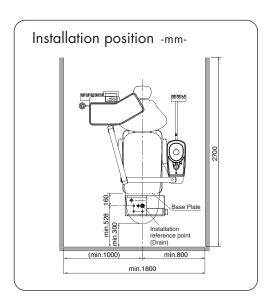

The manufactures reserve the right to alter and improve the specification without notice.









Belmont combines its dental industry knowledge with innovative technology.

efficiency


Through the perfect balance of an ergonomically designed chair and a unit offering smart user interfaces and operational ease, the Clesta ellI creates the ideal dentistry environment for both patient and dentist and ensures efficient and effective practices.

Designed to deliver
unparalleled patient comfort
in any position, the Clesta elll
chair features an axis backrest,
an ultra-smooth tilting mechanism,
and ergonomically shaped
cushions.

Crgonomics


Inheriting the easeenhancing tradition of Clesta dental systems, the Clesta eIII offers a smooth-moving operator's table and easy-access handpiece holders to minimize stress during use.

Case

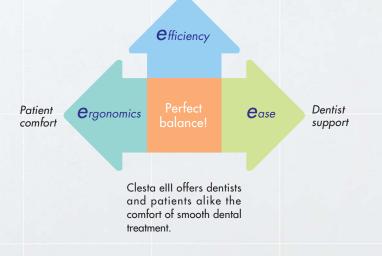





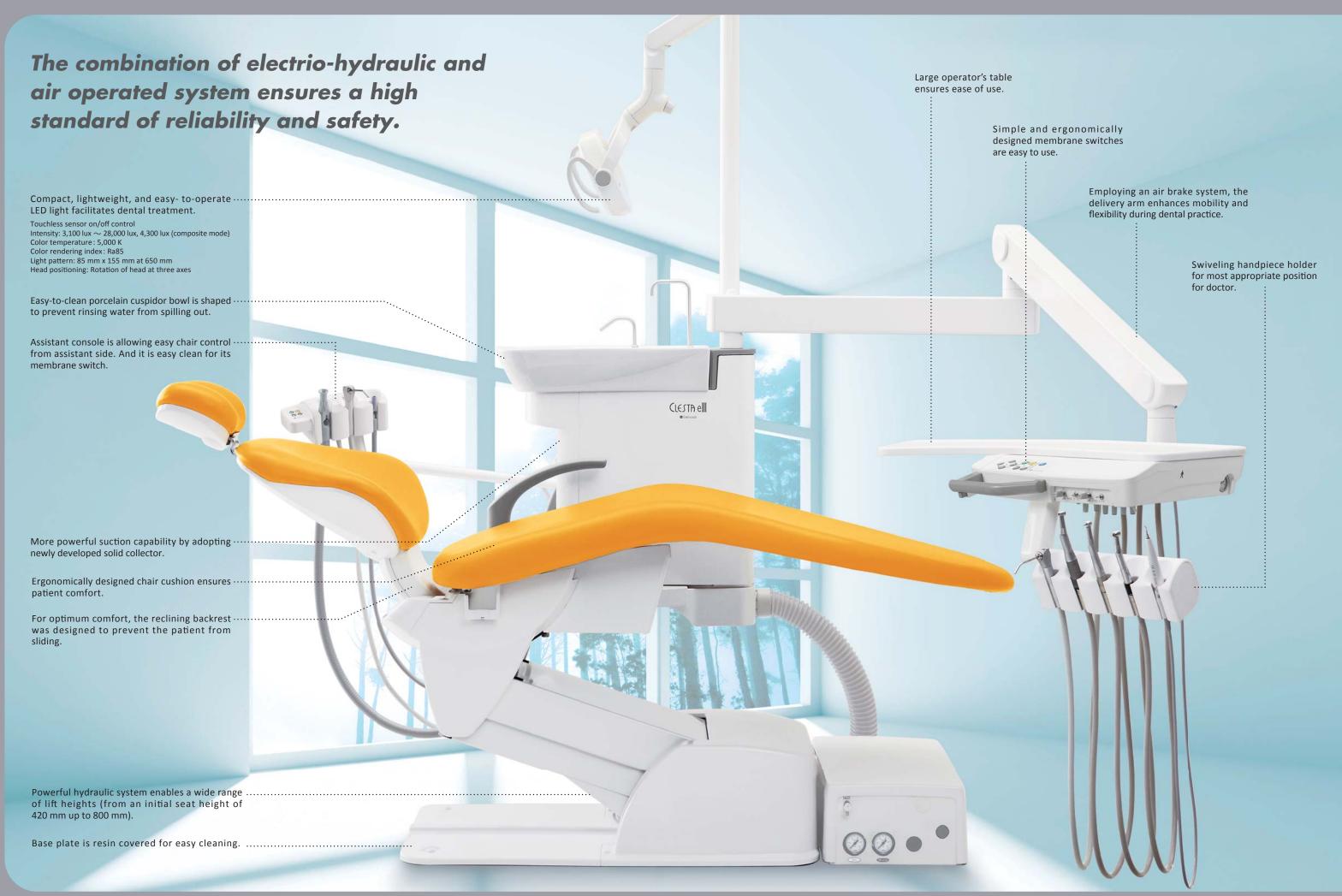

A variety of optional accessories provides greater flexibility of dental treatment.



Flushing system / Flushout tray

The flushing system enables easy cleaning of handpiece water lines. The flushout tray as an extra option allows hands-free, flushing out.



For use with the cupfiller, the water heater helps protect temperature-



Assistant sub-tray increases assistant's work efficiency. (Metal sub-tray is not included) Moreover, curing light holder enhances the operability of the curing light HVE tips and SE tip are optional.



Service outlet / Quick connector

Air and water service outlets and quick connectors for external instruments bring added convenience.



Double articulating headrest

Suitable for children and adults, the double articulating headrest adjusts easily for improved patient



Bracket for flat-screen monitor

This LCD (Liquid Crystal Display) bracket can be adjusted to suit individual needs (monitor not included)

Multifunction foot control

The foot control contains two switches

control switch for up, down, forward,

programmed-position control switch for 0, 1, 2, and last positions.

and backward movement and a



io5000 light

Touchless sensor on/off control Intensity: 28,000 lux (High), 22,000 lux (Medium) & 15.000 lux (Low) Color temperature: 5,000 K Light pattern: 72 mm x 194 mm at 700 mm

900000



048 light

Intensity: 28,000 lux (High), 18,000 lux (Low) Color temperature: 3,700 K Light pattern: 80 mm x 200 mm at 725 mm



Bottled water system

The bottled water system offers a choice of 750-ml or 2-liter bottles and can be integrated into the unit to supply handpiece water. The aesthetic cover is an extra option



Panorama-size film viewer

The Dental-size or panorama-size film viewer features an LED backlight.



Operator table accessories-waste receptacle and cotton container with a stand-support hygienic practices.



Contorol panel for electric scaler and micro-motor.



A swing-out armrest for the right side promotes patient comfort and



Pedodontic chair

Clesta elli PEDO has a smaller backrest, a narrower seat and a cuspidor positioned lower to accommodate young children. The double articulating headrest and the armrest on the right side

Upholstery Color Options

MR1 DARK BLUE

CG18 GRAY TEAL



MR2 PURPLE



CG19 BLUE



LIGHT GREEN



CG22 LIGHT PINK



SAFFRON

CG17



CG15 LIME GREEN



CG2



BEIGE















ES11





