



The sterilization company

Committed to quality and innovation

# Newseal



# Newseal

## The professional sealer



The NEWSEAL system enables surgical instruments and handpieces to be wrapped before being put in the autoclave so as to maintain sterility after sterilisation has been completed.

The new sealing system uses a self-regulating heating element which produces a 10 - 11 mm seal with no risk of accidentally burning the paper.

The 10-11 mm seal operated by the NEWSEAL further accomplish the seal thickness requested by the norm EN-ISO 868-5 (6 mm) and it grants the higher safety standards.


The technology employed by the NEWSEAL system ensures great reliability and easy and safe use even by any member of staff.


The state of the art materials and strict quality control in production and testing provide high performance and the utmost reliability.

- Automatic 10 - 11 mm seal
- Self-regulating heating element
- Scorch protection


NEWSEAL PLUS is equipped with an additional visual indication and a sound signal that will inform the user that sealing have been finished as requested by the norm EN-ISO 11607-2.

NEWSEAL ...the professional choice

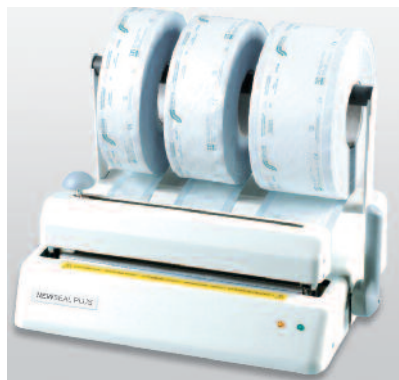

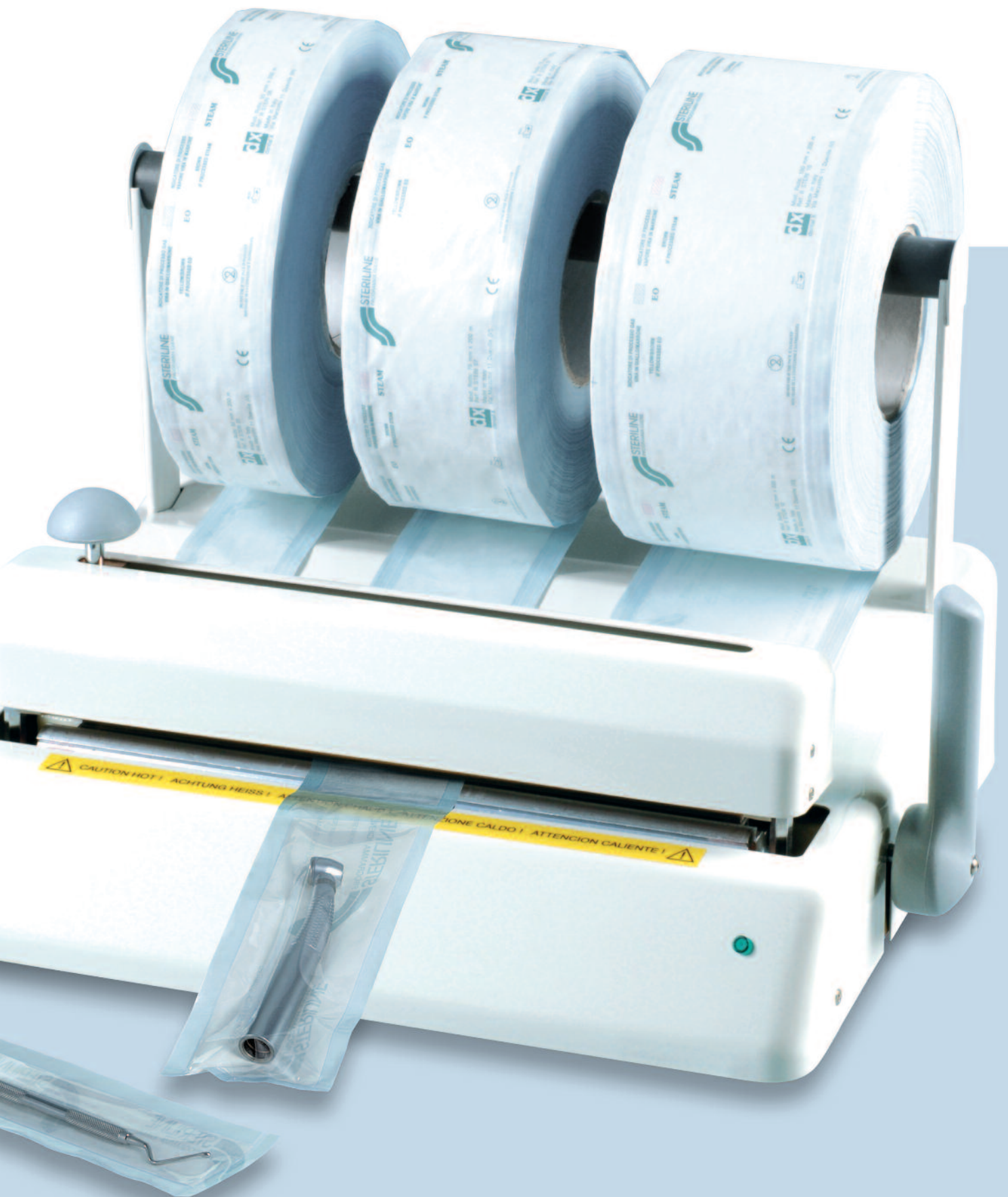

	Newseal	Newseal Plus
Sealing	10-11 mm	10-11 mm
Blade snapback mechanism after cutting	■	■
Automatic start of sealing as soon as the lever is down	■	■
Visual display after seal	/	■
Reel support	■	■
Lever scorch protection	■	■
Metal cover	Bright White 9016	Bright White 9016
Dimensions (HxDxW)	335x323x421	335x323x421
Weight	8 kg	8 kg
Electric Power	110-240 V	110-240 V
	50-60 Hz	50-60 Hz
Consumption Max	400 W	400 W
Wall support	■	■
Pre-heating time	1,5 min	1,5 min
Heating system	Self-adjusting	Self-adjusting
Visual indication	1 Led	2 Led
CE	■	■
ISO 11607-2	/	■
ISO 868-5	/	■

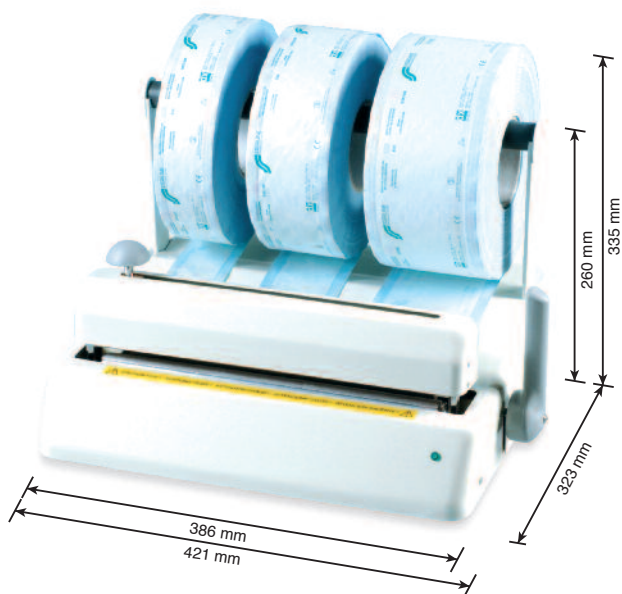

Dental X is a certified company ISO 9001/Vision 2000 and all products are released complete with compliance certificate and statement.



Dental X reserves the right to make alterations without prior notice.

always better



dental X spa  
 via marzotto 11  
 36031 dueville vicenza  
 tel. +39 0444 367400  
 fax +39 0444 367436  
 e mail dentalx@dentalx.it  
 internet www.dentalx.it